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ETB-1010-BW

eldonex[®]

THANK YOU FOR BUYING AN ELDONEX PRODUCT.

Do you want to be the first to know about other news of Eldonex? Follow us on social media.

#eldonex

- Before putting this product into operation please read through the User's Manual carefully, even if you are already familiar with the use of similar products. Use the product only as described in this User's Manual.
- Keep this User's Manual for future reference. It is an integral part of the product and may contain important instructions for commissioning this product, its operation and cleaning.
- Make sure that all other people handling this product are familiar with this Manual. If you hand the product over to other people, make sure that they read this User's Manual, which shall be given to them together with the product.

We recommend keeping the original packaging of the product, proof of purchase and the warranty card, if supplied, at least for the warranty period. In case of transport, we recommend packing the product in the original packaging, in which it has been delivered in, as it protects it best against damage during transport.

Scan this QR code to get the latest version of the User's Manual. In your smartphone, open the QR code reader application and point the smartphone at this code – the page for downloading the latest version of the User's Manual will open:

You can also find the latest version of the User's Manual on our website www.eldonex.com

The text and technical specifications are subject to change.

1 Important Safety Warnings and Instructions

Important safety warnings and instructions

- This appliance is intended for use by children 8 years old and older and persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, provided that they are under supervision or have been given instructions concerning use of the appliance in a safe manner and they understand the potential hazards.
- Children must not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.
- Unpack the product carefully and be careful not to throw away any part of the packaging material until you find all the parts of the product.
- Remove all packaging and marketing materials from the product before first use.
- Take extra care if the product is to be used near children. Always place the product out of their reach. Children should be supervised to ensure that they do not play with the product.
- Do not allow the product to be used as a toy.
- Do not use the product other than as described in this User Manual.
- The product is intended for indoor use; do not use it outdoors.
- Never use accessories that are not supplied with this product or are not intended for this product.
- You should not repair or modify the product yourself under any circumstances! Entrust all repairs and adjustments to this product to a professional company / service company. Interference with the product during the warranty period exposes you to the risk of losing warranty claims.
- Use and store the product away from explosive gases, combustible, flammable and volatile substances (solvents, varnishes, adhesives, etc.).
- Do not expose the product to extreme temperatures, direct sunlight, excessive humidity or excessive dust.
- Do not place the product near heaters, open fire or other appliances or devices that are sources of heat.
- Do not use the product for purposes other than those for which it is intended.
- Do not spray the product with water or any other liquid.
- Do not shake the product and make sure that it does not fall to the ground.
- The manufacturer shall not be liable for damages caused by incorrect use of the appliance and its accessories.
- If you do not use the product for an extended period of time, remove the batteries from the product, otherwise they may leak into the battery compartment and irreparably damage the product.

- Do not clean any part of the product in the dishwasher.
- This appliance has only been designed for cleaning teeth, gums and tongue.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- Stop using the appliance and consult your dentist/doctor if excessive bleeding occurs after use, if bleeding continues to occur after 1 week of use or if you experience discomfort or pain.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use this appliance.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use this appliance.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.

2 Packaging contents

- 2x sonic toothbrush
- 4x standard brush head
- 4x sensitive brush head
- 2x travel case
- 2x wireless charging base
- User manual


3 Technical specifications

- Sonic technology
- Sonic vibration frequency: 40 000 times/min.
- 5 cleaning modes
- 500 mAh Ion-battery
- 2 minutes automatic timer
- 30 seconds intervals
- Waterproof IPX7
- Power supply: 5V ~ 1A (USB-C)
- 40 uses on a single charge

4 Description of each part


A

- 1 Sonic toothbrush
- 2 ON/OFF/modes button
- 3 Modes
- 4 Power indicator light
- 5 Brush head
- 6 Wireless charging base

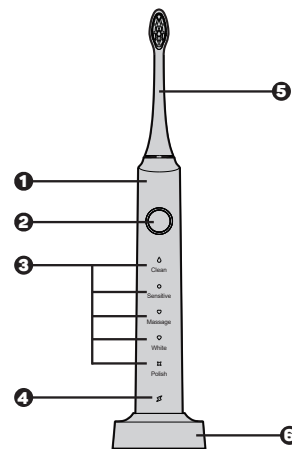5 First steps

- 1 Align the brush head so the bristles point in the same direction as the front of the handle.
- 2 Firmly press the brush head down onto the metal shaft until it will not go any further.

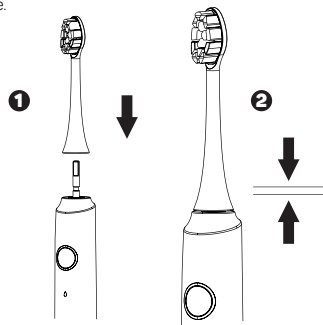NOTE:

It is normal to see a slight gap between the brush head and the handle.


6 Overview

- 1 Wet the brush head and apply toothpaste to it.
- 2 Place the toothbrush on the teeth at a slight angle (45 degrees), pressing firmly to make the bristles reach the gumline and then turn it on.
- 3 Press the button to select one of the five cleaning modes. Press the button for less than 5 seconds or the brush will turn off.
- 4 Move the bristles back and forth slowly and slowly move between teeth. At the same time, gently push the brush against the tooth.
- 5 For maximum efficiency, the cleaning is divided into 4 sections. Each section is given 30 seconds. After a 30 second cycle, the brush pauses and alerts you to move to the next section. After 2 minutes the brush switches itself off.
- 6 Rinse the head and body of the brush under running water and then wipe them with a cloth or a towel.
- 7 Place the body back into the base or store it in a travel case.

NOTE:

The bristles should slightly flare. It is not recommended to scrub as you would with a manual toothbrush..

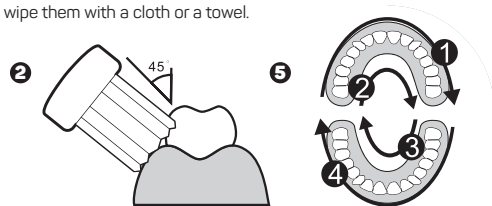7 Modes

- This product has an intelligent memory function, automatically entering the previous usage mode when activated next time.
- Each mode works for 2 minutes.

CLEAN	Basic mode for perfectly cleaned teeth
SENSITIVE	Gentle cleaning for sensitive teeth and gums (also suitable for users not used to sonic cleaning)
MASSAGE	Very gentle stimulation and massage of teeth and gums
WHITEN	Switching between cleaning and polishing modes continuously helps to whiten and maintain the surface of teeth
POLISH	Ultra high frequency of cleaning to polish and brighten teeth

8 Charging

1. Plug the cable into an electrical socket and insert the toothbrush into the charging base. The charging indicator lights up. If you are charging your device for the first time, full charging will take 12 hours.
2. Further charging of the discharged battery to fully charged will take 8 hours. In normal use, a fully charged battery lasts 20 days. The red light on the battery indicator lights up when the battery is running low. Whilst charging, the indicator blinks red and turns green when it is fully charged.
3. If the device is not used for a long time, it is normal that the indicator does not light up immediately. The indicator lights up within 30 minutes.

9 Maintenance

Before cleaning, make sure the appliance is not plugged in.

Toothbrush:

1. Remove the brush head.
2. Rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste.
3. Wipe the entire surface of the handle with a damp cloth. Do not immerse the entire body of the brush in water.

**NOTE:**

Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage. Do not put the brush or its parts in the dishwasher.

Brush head:

1. Rinse the brush head and bristles after each use and air dry naturally.
2. Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week. Rinse the travel cap as often as needed.

Charging base:

1. Unplug the charging base before you clean it.
2. Wipe the dirt from the charging stand with a soft cloth. Do not immerse the charging base in water.

Travel case:

1. Wipe the surface and inside of the travel case with a damp cloth.

**NOTE:**

It is recommended to replace the brush head every three months (once the bristles fades), and replace it when the bristles are bent or deformed.

10 Notices

INSTRUCTIONS AND INFORMATION REGARDING THE DISPOSAL OF USED PACKAGING MATERIALS

Packaging materials and old appliances shall preferably be recycled. The packaging materials of this product can be disposed of as sorted waste. The same applies to the plastic bags made of polyethylene (PE) and other materials – please hand them over for recycling.

DISPOSAL OF WASTE ELECTRICAL & ELECTRONIC EQUIPMENT

In accordance with the European Directive 2012/19/EU on waste electrical and electronic equipment (WEEE), this appliance is labeled with this symbol on the product or its packaging to indicate that this product shall not be treated as household waste. It must be disposed of at a collection site for recycling of waste electrical and electronic equipment. The proper disposal of this product helps to prevent adverse impacts on the environment and human health, which could otherwise result from improper waste disposal. The disposal shall be carried out in accordance with waste management regulations. For more detailed information about the recycling of this product, contact the local authorities, household waste disposal services or the shop where you have bought the product. For proper disposal, renewal and recycling, hand the products over to the designated collection sites. Alternatively, in some EU or other European countries, you can return the products to your local dealer when buying an equivalent new product. The proper disposal of this product helps to preserve valuable natural resources and prevent potential adverse impacts on the environment and human health, which could result from improper waste disposal. For details, contact your local authorities or the nearest collection facility. In case of improper disposal of this type of waste, fines may be imposed in accordance with national legislation.

For business entities in the EU countries

If you want to dispose of electrical and/or electronic equipment, please contact your dealer or supplier for the necessary information.

Disposal in other countries outside the European Union

This symbol is valid in the European Union. If you want to dispose of this product, contact your local authorities or your dealer for the necessary information about the proper way of disposal.

In accordance with the law, we provide the take-back, processing and free environmentally friendly disposal of electrical and electronic equipment through the collective system of REMA Systém, a.s. The main purpose is to protect the environment by ensuring the efficient recycling of waste electrical and electronic equipment.

For the current list of collection sites, see the web www.rema.cloud.

This product meets all the basic requirements of the applicable EU directives. It is CE marked. This mark indicates the conformity of the product's technical characteristics with the applicable technical regulations.

The EU Declaration of Conformity is available at www.eldonex.com.

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Eldonex[®] is a registered trade mark